

Thermostatic water tank



Features


Microprocessor control with timing and temperature correction function.

- The shell is made of high quality steel plate ,surface coated.
- 304 stainless steel chamber, easy to clean and anti rust.
- Between the inner and outer space are filled with high-quality insulation materials,anti-corrosive and good appearance.
- With timing function
- LCD display

Specifications

Model	LTT-600 (3 Usage)	LTT-420 (3 Usage)	LWT-600	LWT-420
Temperature Range (℃)	RT+5~100		RT+5~95	
Temperature Fluctuation (℃)	±0.5			
Temperature Uniformity(℃)	±0.5			
Power Supply	220V / 50Hz			
Power Rating (W)	2000±10%	1500±10%	500±10%	300±10%
Chamber Size (W×D×H)cm	60×30×15	42×18×11	60×30×15	42×18×11
Volume(L)	27	8	27	8
Package Size (W×D×H)cm	81×48×48	65×35×34	81×48×47	65×34×34
Net/Gross Weight (kg)	19/33	11.5/13	19/33	11.5/13

LTF-8D Water Bath

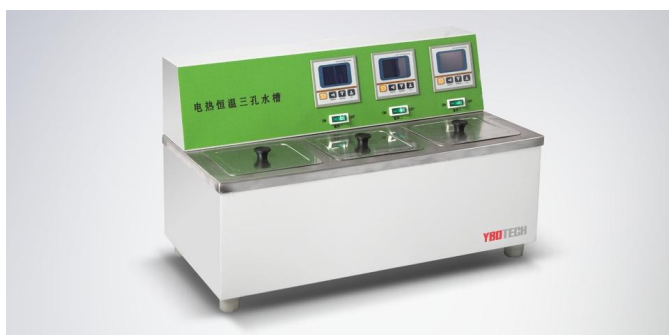

Features

- Stainless steel chamber and top cover.
- 3 independent water tanks and temperature controls.Can separately set and control temperature.
- Freely set operating temperature, with alarm function.
- LCD display for set temperature and measured temperature .
- Equipped with leakage protection.

Specifications

Model	Power Supply	Temperature Range (℃)
LTF-8D	220V,50HZ	RT+5 ~ 99
Power Rating (W)	Temperature Range (℃)	Chamber Size (W×D×H)cm
800	±0.5	(14×16.8×11.5)×3
Volume(L)	Package Size (W×D×H)cm	Net/Gross Weight (kg)
8	58×35×38	9/14